# insta-FRESH® Dental Appliance Cleaner

DIRECTIONS FOR USE

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MODE D'EMPLOI

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
INGREDIENTS: Sodium Bicarbonate, Citric Acid, Sodium Carbonate, Potassium Persulfate Compound, Corn Syrup Solids, Sodium Percarbonate, Sodium Sulfate (may contain silica), Sorbitol, Tetraacetylenediamine (TAED), PEG-180, Sodium Lauryl Sulfoacetate, Flavor, PEG-8, Magnesium Stearate, FD&C Blue #1 Aluminum Lake (CI 42090), FD&C Blue #2 (CI 73015), FD&C Yellow #5 Aluminum Lake and FD&C Yellow #5 (CI 19140)

Manufacturer:  
Dentsply International Raintree Essix  
7290 26th Court East  
Sarasota, FL 34243 USA  
800-866-6884  
[www.smileshop.com](http://www.smileshop.com)

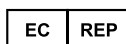Emergo Europe  
Prinsessegracht 20  
2514 AP, The Hague  
The Netherlands

## Insta-FRESH® Dental Appliance Cleaner

**INTENDED USE:** Product is intended to clean a variety of removable dental appliances.

**INDICATIONS FOR USE:** Product is indicated for patients who are using removable dental appliances such as retainers, aligners, dentures, sports mouthguards, nightguards, bruxing and TMJ appliances, which need to be cleaned at home.

**CONTRAINDICATIONS:** This product is contraindicated for patients allergic to Persulfates.

**⚠️ WARNINGS:** Contains persulfates which are a known allergen. Keep out of reach of children. Use only as directed. If product comes into contact with eyes, wash eyes thoroughly with water for several minutes. Do not put tablets or cleaning solution in mouth. If ingested, seek medical assistance immediately. In the USA, please call Poison Control Center at 1-800-222-1222. In the UK, please call NHS Direct at 0845 4647.

**⚠️ PRECAUTIONS:** Do not allow wire appliances to soak for more than 3 minutes. Remove promptly and rinse thoroughly with water. Wash hands thoroughly before and after handling tablets. Do not mix with other cleaners. Single use only. If re-used, product will not perform as intended. Store in a cool dry place.

**⚠️ CAUTIONS:** **Keep out of reach of children. Do not put tablets or cleaning solution in mouth.**

**ADVERSE REACTIONS:** Allergic reactions.

**STEP BY STEP INSTRUCTIONS:** Use daily for best results.

**1** Drop one tablet into enough warm (not hot!) water to cover appliance.

**3** After 3 minutes\* appliance is clean and odor free!

**2** Place appliance into bubbling green solution.

**4** Rinse appliance and hands thoroughly with running water.

Kills 99.99% of common odor-causing bacteria including: Staph aureus, E. coli, P. aeruginosa and F. nucleatum. Data on file. Independent 2015 U.S. lab testing (ASTM method E2315-03).

**\*DO NOT allow wire appliances to soak for more than 3 minutes. Remove promptly and rinse thoroughly with water.**

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